

Warnings:

1. Do not use this device if you have poor blood circulation, diabetes or related conditions.
2. Do not use this device on skin that is inflamed, infected, wounded or sensitive (such as facial skin).
3. Do not use this device if any component becomes damaged, dented or rusted.
4. Consult a physician or skin specialist if you have any pre-existing skin conditions before using this device.
5. If the device causes irritation or other undue stress to your skin, seek medical attention immediately.
6. Keep the device away from young children. It is **NOT** intended for use by children.
7. Do not use or store the device near water or liquid, and ensure that the targeted area of skin is dry and clean before use.

Specifications:

Power Input	AC/DC 100-240V, 50/60Hz, 0.2A
Power Output	12V DC Motor, 1.0A
Dimensions	95 x 70 x 70mm
Weight (gross)	18 oz
Weight (net)	14 oz (including Power Adapter)

Warranty:

The manufacturer warrants the product to be free from defects in workmanship or material under normal use for six (6) months from the date of purchase. The manufacturer's obligation under this warranty is limited to replacing, adjusting or repairing the unit if returned along with proof of purchase. This warranty is void if the unit has been tampered with or maliciously damaged.



CallusClear USER'S MANUAL

The CallusClear System

Callus Removal Pedicure Solution

www.HomeCareLogic.com


Technical Support: Distributed by HomeCareLogic, LLC
21 Grand Avenue, Suite 124, Palisades Park, NJ 07650
Toll Free: 1-877-606-8080
<http://www.HomeCareLogic.com>
support@homecarelogic.com



CallusClear Callus Removal Pedicure Solution

The CallusClear is the most effective callus removal/pedicure solution available today.

The CallusClear can eliminate even the worst callus problems safely and quickly, or give you an easy and convenient way to cure rough skin patches for beautiful, healthy skin.

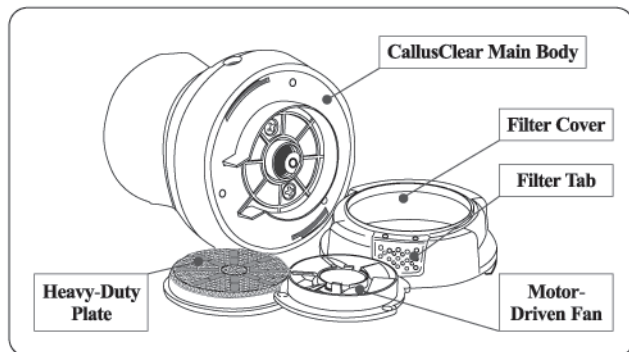

Operating Instructions:

1. Ensure the switch on the side of the CallusClear is in the Off position.
2. ***Remove the outer lid from the device by pulling straight up; if you twist when lifting the lid, the filter cover may become detached.**
3. Ensure that the grinding plate and filter cover are securely in place.
4. Adjust the switch to your desired spin strength setting (Lo or Hi).
5. Grip the device firmly and apply the spinning grinding plate to the callus (or rough skin surface). The device should be continuously rotated in small clockwise circles while the plate is against the targeted surface.
***Do not press the device continuously in a fixed position against one area of skin, as this may cause a buildup of heat and irritate that area.**
6. When you are finished using the device, turn the power off and unplug the cable from the power outlet. You can open the filter tab alongside the spinning plate and tap the collected residue into a wastebasket.

(Optional) Dismantled Cleaning Instructions:

1. ***Remove the outer lid from the device by pulling straight up; if you twist when lifting the lid, the filter cover may become detached.**
2. Remove the filter cover (that is surrounding the spinning plate) by gently twisting it counterclockwise and lifting up.
3. Remove the spinning plate attachment; there are two indented portions around the edges of the fan underneath where you can grab the edges of the plate attachment to lift up easily. (When re-attaching the plate, make sure that the three orientation nubs are aligned with the three holes along the outer edge of the fan.)
4. The loose portions (filter cover, plate and fan attachments) can be cleaned using any method that does not interact with plastic or aluminum, and are dishwasher-safe. The main body of the device should only be cleaned with a dry cloth and should not be exposed to any liquids or chemicals.